


# LEVITY

by ORA Medical, Inc.

Download Ora Connect for:

- Levity video guides
- Suggested workouts
- Progress tracking

Available on iOS & Android!

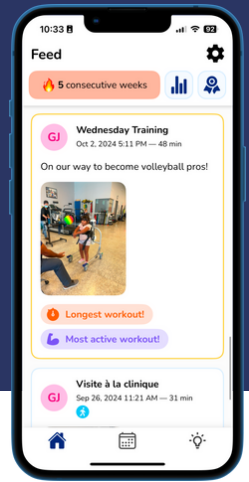

## Unpacking the Levity



- Remove the device.
- Check your order.

Tip: Lay the package on its side to pull the device out with ease



- Loosen the levers.
- Unfold the Levity.

Tip: Hold the bottom frame down then pull the upper frame up.




- Press down on the sliders to ensure engagement of safety notch.
- Tighten both the levers.



- Slide the hip cushions onto the harness rail.
- Tighten the knobs in position.

## Changing Weight Modules

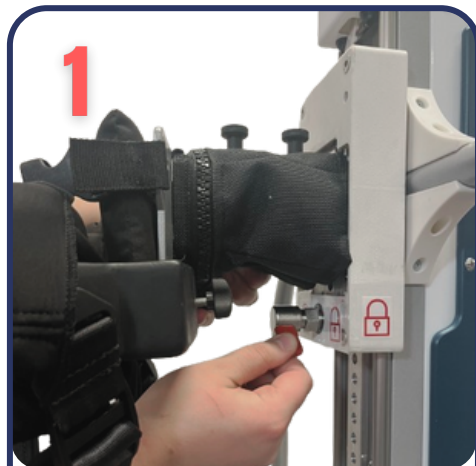

- Lift the harness to the top.
- Lock the brake in the highest position, by rotating the lever.



- Lift the rear cover
- Detach the buckle.



- Slide the housing out.
- Slide the new housing in. A “click” should be heard.
- Securely attach the buckle and close the rear cover.



Refer to the user manual.

## Contact Us


+1 (514) 839-1773  
info@oramedical.ca  
80 rue Queen, Montréal, QC  
H3C 2N5, Canada

## Thank you for choosing our product !

Our team is excited to accompany you as your child takes their first steps with Levity!

Share your experience with us on Instagram by tagging @ora\_medical!

## Installing the User

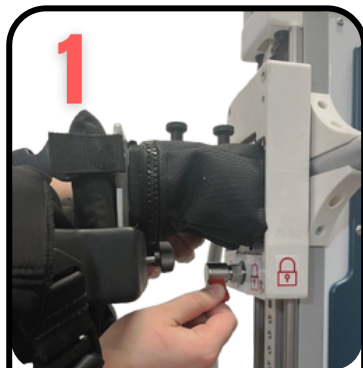

- Adjust the harness to hip level.
- Lock the brake by rotating the lever.
- Unclip the front buckles on the pelvic support



- Place the user in the device.
- Pass the straps between their legs.
- Attach the front buckles.

*Tip: The belt can be installed on the user prior to installation in the device*



- Adjust front and rear straps to enhance comfort and posture of the user.

*Tip: Tightening the straps prior to installing the user, makes adjustments easier.*



- Buckle the hip belt around the user's waist.
- Adjust the straps for improved comfort.

## User Adaptation



### Swivel Control

High Position: Full swivel  
Mid Position: Partial swivel  
Low Position: No swivel



### Speed Control

- Turn the handle clockwise to limit speed.

*Tip: Use the same adjustment on both wheels to avoid affecting the direction of movement.*



### Torso Rotation Control

High Position: Allows torso rotation.  
Low Position: Blocks torso rotation.



### Side Inclination Control

High Position: Allows side inclination.  
Low Position: Blocks side inclination.



Refer to the user manual.